




J-TALK™

TALKING BATHROOM SCALE

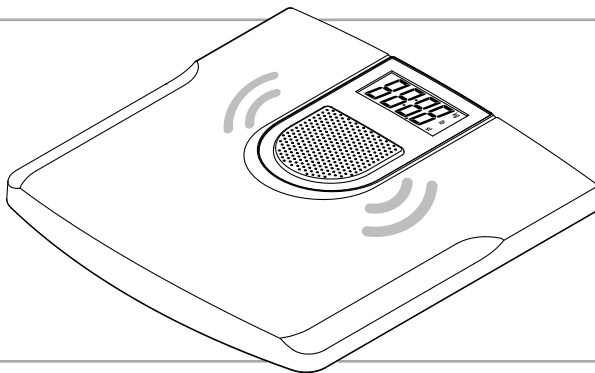

Thank you for purchasing the J-Scale® J-Talk™ digital scale. Please read all operating instructions carefully before use. This electronic scale is a precision instrument. With normal care and proper treatment, it will provide years of reliable service.

For further information or troubleshooting, please visit our website at www.jscale.com

FEATURES & SPECIFICATIONS

Capacity	200kg / 440lb / 31:7st
Weight graduation	0.1g / 0.2lb
Units	kg / lb / st:lb
Temperature Range	0°C - 40°C
Power source	1 x 9V battery
LCD/Talking	Simultaneous
Auto Off	8 seconds
Volume Control	High, Middle, Low
Low Voltage Message	<i>L</i>
Error/Overload Message	<i>Err</i>



OPERATION INSTRUCTIONS

Weighing


- 1) Place the scale on a flat hard floor.
- 2) Use your foot to tap the scale platform once and turn the scale on.
- 3) A voice message will tell you to wait.
- 4) When the display shows "0.0", a voice message will tell you to step on the scale.
- 5) Now step on the scale and stand still.
- 6) Your weight will be displayed on the LCD and a voice will tell you the weight at the same time.

Caution


To avoid water from getting into the holes and damaging the speaker, please don't step on the scale when you are wet.

Changing the battery


Ensure +/- terminals are correctly connected.

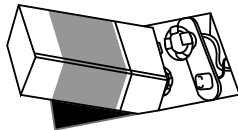

DESCRIPTIONS

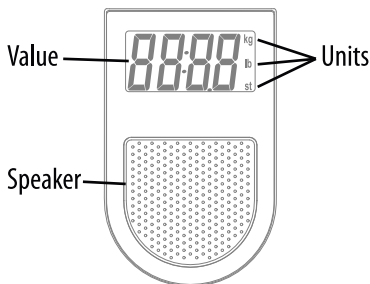

Volume Control
(bottom of scale)

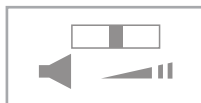

Unit Selector
(bottom of scale)



NON SLIP STRIPS

Included with the scale are non-slip strips to be installed onto the scale's platform. These help prevent your feet from slipping off of the weighing surface, and possible slip-and-fall accidents. Note that we strongly urge you to never step on the scale with wet feet regardless of use of these strips, to prevent accidental slippage or an accident.

DANGER - SCALE SURFACE MAY BE SLIPPERY WHEN WET