


ProScale®

MOUSE

SAFE & SCALE

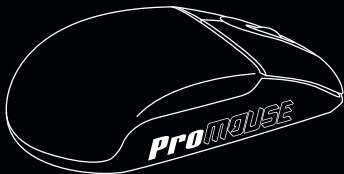

USER MANUAL

Thank you for purchasing the PRO-MOUSE™ scale. Please read **all** operating instructions carefully before use. This **electronic scale** is a precision instrument. With **normal care** and proper treatment, it will provide years of **reliable service**.

Never **load the scale** with more than the maximal capacity. **Overloading will permanently damage it!** Avoiding any exposure to extreme heat or cold, your scale works better when operated at normal room temperature. **Always allow the scale** to acclimate to normal room temperature for at least one hour before use. Give your scale sufficient warm up time. Usually 30–60 seconds before calibration to give the internal components a chance to stabilize. Keep your scale in a clean environment. Dust, dirt, moisture, vibration, air currents and a close distance to other electronic equipment can **all** cause an adverse effect on the **reliability and accuracy** of your scale. **Handle with care**. Gently apply **all** items to be weighed onto tray top. Although this scale is designed to be quite durable, try to avoid rough treatment as this may permanently damage the internal sensor. Avoid shaking, dropping or otherwise shocking the scale. This is a precision instrument and must be handled with extreme care.

OPERATION

1. Place scale on a horizontal flat surface, press [ON/OFF] key to turn on the scale.
2. Wait until "0.0" or "0.00" is displayed
3. Put the object(s) on the weighing platform
4. Using the [M] key, you can switch among the weighing units.

TARE

1. Turn on scale as described above
2. Place the "tare item" on the platform
3. Press [TARE] and wait until "0.0" is displayed
4. Add the "net-weight-item"

CALIBRATION

1. Turn on the scale by pressing [ON/OFF] key.
2. Wait until the LCD displays "0.00" and then press and HOLD [M] key until the "CAL" will be displayed, and then release the [M] key.
3. Press [M] key again, the LCD will display "CAL" followed by the required calibration weight.
4. Add the required calibration weight, the display will show "PASS" and return to normal weighing mode.
5. Put away the calibration weight and press [ON/OFF] to turn off the scale and the calibration is finished!

**PLEASE CHECK THE AVAILABILITY OF WEIGHTS
AT THE STORE YOU BOUGHT YOUR SCALE FROM.**

SPECIFICATIONS

Tare Range	Full capacity
Auto Off	60 seconds
Op. Temp	10-30 °C
Power	2x AAA
Scale Dim.	59x59x20mm
Tray Dim.	57x39mm
Mouse Dim.	120 x 64 x 35mm

MOUSE 100

Capacity	100g
Readability	0.01g
Modes	g, oz, gn, ct
Calibration	100g

MOUSE 500

Capacity	500g
Readability	0.1g
Modes	g, oz, gn, ct
Calibration	500g



ProScale® **MOUSE**

Manufactured in China - Fabriqué en Chine